cálido confort

Bota duoterm  
Duoterm boot  
Bota duoterm

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**WARNING**

**Please read the instructions for use carefully prior to using the product. Store these in a safe place for future reference.**

The Duoterm boot contains a removable heating pad from the DAGA PVC line in its interior with a plastic finish that can be used separately to the boot itself.

The Duoterm boot can only be used with an EP model pad in its interior.

Prior to connecting the pad, check that your power supply voltage is identical to the one indicated on the label included on the product or on the plastic print embedded onto the surface of the pad body. The mains connection cable must not be tangled or wrapped around the pad during use.

Do not pull on the cable to unplug it or use it to carry the pad.

Do not wrap the cable around the control.

Do not cover or lay the regulator control on top of the heating pad during use in order to ensure its proper operation.

To ensure the proper use and operation of the Duoterm boot, do not cover it with large or heavy objects or lay on top of other heating elements.

Ensure it is not folded during use.

Do not stick any pins or needles or any other pointy metal objects into the Duoterm boot.

This product has been specifically designed to be used directly on the floor. However, regular footwear should not be worn while using this device.

Do not use on a wet or damp floor.

Ensure that the Duoterm boot does not get wet.

Never immerse the boot in water or use close to the same.

If any parts of the Duoterm boot get wet or damp, do not use the boot until it is dry.

Refer to the maintenance section of this manual for instructions on how to clean the device.

People that are disabled, handicapped or that cannot feel heat, as well as children or sick people should only use this product under direct adult supervision.

Children should be supervised at all times while using this product.

This product is not a toy.

It should not be used under the effects of alcohol, sleeping tablets or sedatives.

## ENGLISH

Do not use the Duoterm boot to heat up animals.

While the electromagnetic field generated by its heating pad is extremely low, we recommend checking with your medical specialist prior to using this product if you are fitted with a pacemaker or any other similar cardiovascular device.

While the product has a temperature reduction and automatic disconnection function, we still recommend adjusting the temperature controls to the minimum operation level (level 1) during long periods of use.

Do not sleep while the Duoterm boot is connected.

This device is intended solely for domestic use and is not intended for medical use in hospital.

Materials will suffer wear and tear over time and after prolonged use.

Carry out a regular inspection of the Duoterm boot, including the control and flexible cable to check for any signs of wear or damage. Do not use the product if you spot any visible flaws and contact any B&B TRENDS SL Authorised Technical support service.

Unplug the product immediately from the mains in the event of any breakdown and contact an authorised Technical support service. In order to prevent any risk of danger, do not open the control or connection. Only qualified technical personnel from the brand's Official Technical support service may carry out repairs or procedures on the device.

Any incorrect use or improper handling of the Duoterm boot will render the warranty null and void.


Only an authorised technical support service centre may carry out repairs on this product.

**B&B TRENDS SL.** accepts no liability for any damage that may be caused to people, animals or objects owing to the failure to observe this warning correctly.

## ERGONOMIC LINE - FEET CARE

The reference of the Ergonomic Line that you have purchased is found on its original packaging and the textile label attached to the product itself.

### DUOTERM BOOT with EP heating pad

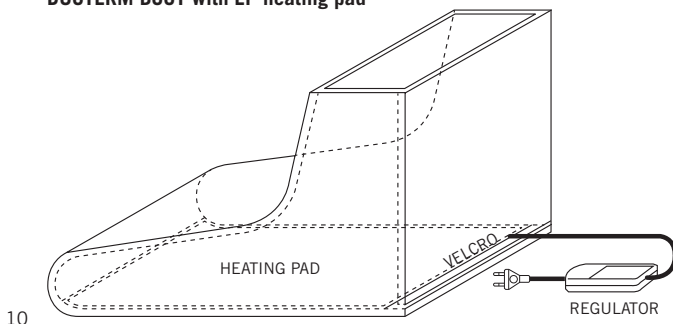10


The anatomic design that has been especially designed for the lower limbs allows the Duoterm boot to mould perfectly around the same.

The application of dry heat generated by the Duoterm boot provides a pleasant sensation of comfort. This product is recommended for people that spend a lot of time sitting down or those that are very sensitive to cold.

The Duoterm boot also contains an EP plastic pad that is removed from the rear pocket of the boot, providing additional benefits of the heating pad separate to the boot itself.

## CONTROL

Control screen for EP model included in the Duoterm Boot:
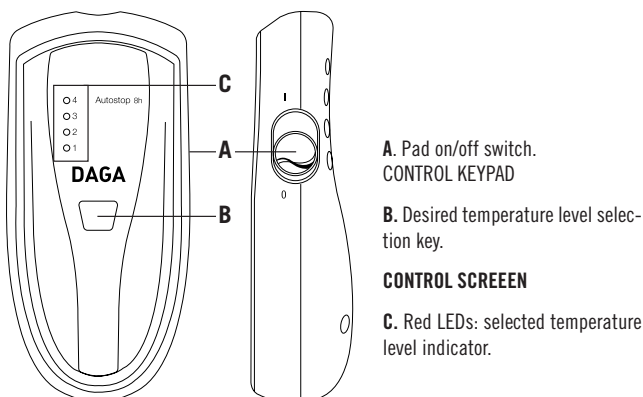
## BASIC FUNCTIONS

The Duoterm boot includes a new generation of controls from the DAGA-TRADITIONAL innovative line with a range of additional functions such as:

- **Temperature selector:** 4 levels, 1 being the lowest and 4 the highest temperature level.
- **Automatic temperature reduction:** This function allows the control to regulate the removable pad use time and temperature parameters for their gradual and automatic control. This enables the pad to operate uninterrupted for 90 minutes at high temperatures (position 3 and 4) after which time, the temperature is automatically reduced, with level 1 lighting up on the screen.

If you wish to deactivate this function, press the central key to allow you to once again select the higher temperature.

- **Automatic safety stop or Autostop:** In order to minimise the risk posed to the user from inadvertently leaving the Duoterm Boot heating pad on for a long period of time, this series includes a special safety automatic stop feature, which activates after the device has been left

## ENGLISH

on for 8 continuous hours if the control key has not been activated. The automatic stop time will go back to the start every time the central key is pressed.

After activating this Autostop feature, the removable pad will remain disconnected and inoperative with the control screen indicates the reason with a flashing red LED displaying the previous temperature level used.

If you wish to deactivate this function, press the central key to allow you to once again select the desired temperature.

### INSTRUCTIONS FOR USE

Prior to connecting the Duoterm boot to the mains, check your power supply voltage to ensure it is the same as the one indicated on the label included on the product.

Connect the product to the proper mains network.

Select the ON position indicated with the “I” symbol located on the side sliding button on the control in order to turn on the Duoterm boot.

Wait until the automatic start check is completed which will finish when the level 1 temperature lights up.

The boot is now ready for operation.

Press the central key several times selecting one of the 4 optional levels in order to set the desired temperature level.

The regulator operates continuously at the temperature level that has been pre-set with the LED lit up on the regulator control indicating the selected temperature.

The maximum operation time is 8 hours, after which an automatic safety stop is activated, displaying the last temperature used flashing on the screen. The Duoterm boot will stop producing heat when in this state.

If you wish to deactivate this function, press the central key to allow you to once again select the desired temperature.

The pad that comes with the boot includes rapid heating technology to allow the heat to be felt only 5 minutes after the control has been started-up, regardless of the temperature selected.

We recommend adjusting the control to the minimum level (level 1) when using the product for prolonged periods of time. However, the pad that comes with the product will automatically reduce in temperature if this is not observed as a safety measure 90 minutes after continuous operation at high levels (positions 3 and 4) to prevent excessive heat exposure.

Select the OFF position indicated with the “O” symbol on the side sliding button on the control in order to turn off the Duoterm boot pad or prior to disconnecting the device.

Let the Duoterm boot pad cool down prior to folding it for storage and do not place objects on top of it.

## TROUBLESHOOTING GUIDE:

### Possible operational problems

PROBLEM	RECOMMENDATION
<p>The regulator control screen does not turn on when the On/Off side switch is activated.</p>	<p>Ensure that there is voltage in the mains being used and that the Duoterm boot is correctly plugged into the same. Check that the On/Off side switch is set to the "I" position. If the problem continues, contact an authorised technical support service centre.</p>
<p>The Duoterm boot removable pad does not heat up.</p>	<p>The safety Autostop may be activated. Check that the boot pad is in operation, checking to see if any of the LEDs on the screen are lit up. If in doubt, press the central key to select the desired temperature again or alternatively, turn the side sliding button off and on again.</p>
<p>After performing the initial check, the control screen is not correctly configured (the position 1 LE does not light up) and the Duoterm boot does not heat up.</p>	<p>The control could be experiencing a problem with start-up, blocking the heating pad in the boot from operating for safety reasons. Re-start the device by turning the control off and on again. If the problem continues, contact an authorised technical support service centre.</p>
<p>After turning it on, the control screen is correctly configured (the position 1 LED lighting up by default) but the removable pad in the boot does not heat up or only heats up a little.</p>	<p>If the removable pad does not heat up, cover it with an item of clothing, set the regulator to the position 4 and wait for 45 minutes. If it still has not heated up after this time, contact our authorised technical support service centre. Checks on the pad cannot be reliably carried out in open spaces or outdoors as this may lead to flawed heat readings.</p>
<p>The control screen shows incomprehensible information. The LEDs will not stop flashing.</p>	<p>The control could have a problem, blocking the heating pad from operating for safety reasons. Turn the control off and on again. If the problem continues, contact an authorised technical support service centre.</p>

# ENGLISH

**There appears to be a reduction in the boot temperature upon completion of the operation time.**

Once the boot is operating for approximately 2 hours at high temperatures, the built-in safety system in the control automatically reduces the heat to level 1.

The Duoterm boot reaches temperatures of between 65° and 75 °C in its interior.

This information is meant as a guide only. The measurements are made under laboratory conditions, according to the procedures described in Standard UNE-EN 60335-2-17.

## PRODUCT DISPOSAL AND MAINTENANCE

### PRODUCT MAINTENANCE

Unplug the Duoterm boot when you do not intend to use it for an extended period of time to prolong its shelf-life.

Store it in its original packaging in a dry place, trying to fold it as little as possible and preventing the formation of creases.

Do not cover the pad with other heated elements during operation. Let it cool down prior to storage.

Do not place any objects on top of the device during storage in order to prevent the formation of creases.

### PRODUCT CLEANING

Do not dry clean or wash the outer part of the Duoterm boot with water.

Do not soak or treat with chemical agents in order to prevent flaws in the product.

Use dry foam to clean the outer textile part of the boot, spreading it uniformly over the product fabric, and then combing it clean length-wise.

Use a damp cloth to clean the plastic pad that comes with the Duoterm boot.

If the boot or any of its parts are amp, let them dry prior to re-connecting it to the mains.

Never iron the pad.

### CLEANING THE COVER.

The Duoterm boot contains 2 washable removable covers: the acrylic fur cover inside the boot and the textile cover that can be used with the plastic pad separately for other uses.

To clean either of the covers, first remove them from the boot and the pad.

Both of them can be machine or hand-washed, using liquid detergent at maximum temperature of 30 °C, selecting the delicate wash program.

We recommend using a wash bag to prevent wear and tear of the fur cover and zipper from catching on other clothes.

ENGLISH

**PRODUCT DISPOSAL.**

Do not dispose of the Duoterm boot at the end of its shelf-life. Contact the organic collection centre closest to your home in accordance with the waste electrical and electronic equipment (WEEE) Directive 2002/96/EC to contribute to the protection of the environment.

Request detailed information regarding their current methods of disposing of packaging and scrappage of used devices from your distributor.

We hope that you will be satisfied with this product.

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**Carefully read the instructions for use.**

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**Do not stick any pins or needles or any other pointy metal objects into the pad.**

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**Children under the age of three are not to use this appliance due to their inability to react to overheating.**

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Lavar a máquina, máximo 30 °C.

**Machine wash at a maximum temperature of 30 °C.**

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**Do not dry clean (dry cleaners or launderette)**

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**This product complies with European Directive 2002/96/EC on electrical and electronic devices, known as WEEE (Waste Electrical and Electronic Equipment), which provides the legal framework applicable in the European Union for the disposal and reuse of waste electronic and electrical devices. Do not dispose of this product in the bin, instead going to the electrical and electronic waste collection centre closest to your home.**

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## WARRANTY REPORT

**B&B TRENDS, SL.** guarantees compliance of this product for the purpose for which it is intended for a period of three years for the entire Spanish territory. In the case of breakdown during the term of this warranty, you are entitled to repair or to the replacement of the product at no charge if the former is unavailable, unless one of these options proves impossible or disproportionate. In this latter case, you can then choose a reduction in price or cancellation of the sale, which must be dealt with directly with the sales vendor. This also covers replacement of spare parts provided that the product has been used according to the recommendations specified in this manual for both cases, and has not been tampered with by any third party that is unaccredited by B&B TRENDS, SL. This warranty does not affect your rights as a consumer in accordance with the provisions of Law 23/2003 dated 10 July 2003 on warranties on the sale of consumer goods in the Spanish territory.

### USE OF WARRANTY

Consumers should contact a technical service centre authorised by **B&B TRENDS, S.L.** for repair of the product as any tampering thereof by persons not authorised by Daga, as well as misuse or inappropriate use, will lead to the cancellation of this warranty.

The warranty must be fully completed and delivered along with the receipt or delivery docket for the effective exercise of rights under this warranty.

Claims must take place within three years from deliver of the project according to the provisions of the warranty, and the consumer must inform the sales vendor of the lack of conformity, provided that such claim is made within the statutory period of warranty

This warranty should be retained by the user as well as the invoice, receipt or the delivery docket to facilitate the exercise of these rights.

For technical service and after-sales care outside the Spanish territory, please submit your query to the point of sale where you purchased the item or by contacting B&B TRENDS, SL. directly.

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